

Nucleocapsid CoV-1/2 sdAb

Cat.No. N3605-DBCO; Single Domain camelid antibody, 250 µg single domain antibody

Data Sheet

| | |
|----------------------------|---|
| Reconstitution/ Storage | 250 µg purified antibody, lyophilized from PBS, fluorescence-labeled with DBCO. Albumin was added for stabilization. For reconstitution refer to the NanoTag reconstitution and storage instructions . Reconstitute immediately upon receipt! Avoid bright light when working with the antibody to minimize photo bleaching of the fluorescent dye. For detailed information, see back of the data sheet. |
| Applications | ICC: yes |
| Label | DBCO |
| Clone | 2A7 |
| Subtype | single domain |
| Immunogen | Recombinant protein corresponding to AA 1 to 419 from SARS-Cov-2 Nucleocapsid protein (UniProt Id: P0DTC9) |
| Specificity | Specific for nucleocapsid from SARS-CoV1 and SARS-CoV2 |
| Matching control | P0101 |
| Remarks | This single domain antibody carries a C-terminal Biotin. |

TO BE USED IN VITRO / FOR RESEARCH ONLY
NOT TOXIC, NOT HAZARDOUS, NOT INFECTIOUS, NOT CONTAGIOUS

Background

The severe acute respiratory syndrome coronavirs type 2 (SARS-CoV2) is an enveloped positive-sense single stranded RNA virus that has been identified in the beginning of 2020 (1). It infects human host cells by docking via its spike protein (S) to the ACE2 surface receptor (2) and can cause mild to very severe and even deadly Covid-19 courses (3).

The very abundant Nucleocapsid or N-protein packages the viral RNA (4) and shares only little homology to other abundant members of the coronavirus family like NL63, 229E, HKU1 or OC43. This characteristic makes it a suitable target to discriminate between Covid-19 and other Coronavirus infections.

Unlabeled variants and several modifications of sdAbs like biotin, fluorophore or DBCO conjugation are available.

In **FluoTag®-Q** each fluorophore is coupled to exactly one FluoTag, which in turn binds to its target molecule in a monovalent fashion. The high binding affinity and a coupling efficiency of > 95% guarantees a highly linear relation between the number of target molecules and the intensity of fluorescence. This enables a direct count of the target molecule of interest. The fluorophore is located exceptionally close to the recognized epitope (< 1.5 nm), which is ideal for all microscopy techniques.

Selected General References

The species Severe acute respiratory syndrome-related coronavirus: classifying 2019-nCoV and naming it SARS-CoV-2. et al. Nat Microbiol (2020) PubMed:32123347


SARS-CoV-2 Cell Entry Depends on ACE2 and TMPRSS2 and Is Blocked by a Clinically Proven Protease Inhibitor. Hoffmann M et al. Cell (2020) PubMed:32142651

Epidemiological and clinical characteristics of 99 cases of 2019 novel coronavirus pneumonia in Wuhan, China: a descriptive study.

Chen N et al. Lancet (2020) PubMed:32007143

The SARS coronavirus nucleocapsid protein--forms and functions. Chang CK et al. Antiviral Res (2014) PubMed:24418573

Access the online factsheet including applicable protocols at <https://sysy.com/product/N3605-DBCO> or scan the QR-code.

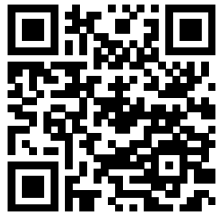

FAQ - How should I store my antibody?

Shipping Conditions

- All SYSY antibodies and control proteins/peptides are shipped lyophilized (vacuum freeze-dried). In this form, they remain stable without loss of quality at ambient temperatures for several weeks.

Storage of Sealed Vials after Delivery

- **Unlabeled** and **biotin-labeled antibodies** and **control proteins** should be stored at **4°C** before reconstitution. **Do not freeze lyophilized antibodies.** Temperatures below 0°C may impair performance.
- **Fluorescence-labeled antibodies** should be reconstituted immediately upon receipt. Long-term storage of lyophilized fluorophore-conjugates may cause aggregation.
- **Control peptides** should be stored at -20°C before reconstitution.

Long Term Storage after Reconstitution (General Considerations)

- **Do not use frost-free (“no-frost”) freezers.** These units periodically warm to remove ice buildup, causing freeze–thaw cycles that can damage antibodies.
- Store vials in areas with minimal temperature fluctuation - preferably toward the back of the freezer, not on the door.
- Aliquot reconstituted antibodies and store at -20°C to -80°C.
- Avoid very small aliquots (<20 µL), as evaporation and adsorption to tube surfaces can reduce antibody concentration and activity.
- Use the smallest practical storage vial to minimize surface area.
- Adding glycerol to a final concentration of 50% prevents freezing at -20°C, allowing storage in liquid form and effectively avoiding freeze–thaw cycles.

Product Specific Hints for Storage

Control proteins / peptides

- Store at -20°C to -80°C

Monoclonal Antibodies

- **Ascites and hybridoma supernatant:** Store at -20°C to -80°C. Prolonged storage at 4°C is not recommended, as proteases present in ascites may degrade antibodies.
- **Purified IgG:** Store at -20°C to -80°C. Adding a carrier protein (e.g., BSA) enhances long-term stability. Many SYSY antibodies already contain carrier proteins - refer to the respective datasheet for details.

Polyclonal Antibodies

- **Crude antisera:** Can be stored at 4°C with antimicrobials added, but -20°C to -80°C is preferred
- **Affinity-purified antibodies:** Less stable than antisera; store at -20°C to -80°C. Adding a carrier protein such as BSA improves long-term stability. Most SYSY antibodies already contain carrier proteins - refer to the respective datasheet for details.

Fluorescence-labeled Antibodies

- Store as a liquid with 1:1 (v/v) glycerol at -20°C, and protect from light exposure

Avoid repeated freeze-thaw cycles for all antibodies!

FAQ - How should I reconstitute my antibody?

Reconstitution

- All purified SYSY antibodies are lyophilized from PBS. To reconstitute the antibody in PBS, add the volume of deionized water specified in the corresponding datasheet. If a larger final volume is desired, first add the recommended amount of water, then adjust with PBS and, if needed, add a stabilizing carrier protein (e.g., BSA) to a final concentration of 2%. Some SYSY antibodies already contain albumin; please take this into account before adding additional carrier protein.

For complete reconstitution, carefully remove the vial cap. After adding water, briefly vortex the solution. To collect the liquid at the bottom of the vial, place the vial inside a 50 ml centrifuge tube padded with paper and centrifuge briefly.

- If desired, small amounts of azide or thimerosal may be added to prevent microbial growth. This is particularly recommended when storing an aliquot at 4°C.
- After reconstitution of fluorescence-labeled antibodies, add glycerol 1:1 (v/v) to achieve a final concentration of 50%. This prevents freezing at -20°C and keeps the antibody in liquid form, effectively avoiding freeze–thaw cycles.
- Glycerol may also be added to unlabeled primary antibodies as a general measure to prevent freeze–thaw damage.
- For further guidance, please refer to our **storage tips** and recommendations for reconstituted antibodies, control peptides, and control proteins.